<u>Global Seafood Alliance Logo</u>

- GOAL Events
- <u>Advocate Magazine</u>
- <u>Aquademia Podcast</u>
- <u>Blog</u>
- <u>Contact</u>
- 0
- f
- 🗙
- in
- 🕨
- <u>Log In</u>
- $\square$ 
  - <u>About</u>
    - Who We Are
    - <u>Our History</u>
    - <u>Our Team</u>
    - Sustainable Development Goals
    - <u>Careers</u>
  - <u>Membership</u>
    - <u>Overview</u>
    - <u>Our Members</u>
    - Corporate Membership
  - <u>Resources</u>
  - <u>Certification</u>
    - Best Aquaculture Practices
    - Best Seafood Practices

Search...

Q

<u>Log In</u>

- <u>About</u>
  - Who We Are
  - <u>Our History</u>
  - Our Team
  - Sustainable Development Goals
  - <u>Careers</u>
- <u>Membership</u>
  - <u>Overview</u>
  - Our Members
  - Corporate Membership
- <u>Resources</u>
- <u>Certification</u>
  - <u>Best Aquaculture Practices</u>
  - Best Seafood Practices
- GOAL Events
- <u>Advocate Magazine</u>

- <u>Aquademia Podcast</u>
- <u>Blog</u>
- Contact



Health & Welfare Health & Welfare

# Molecular biology tools applied to selection programs

Responsible Seafood Advocate logo 1 December 1999 Sam Stern Leonardo Galli

# Searching for genetic markers

Selective breeding is a process which strives to select organisms bearing genes for preferred traits such as rapid growth, disease resistance and fecundity. The decision of which offspring to select for the next generation may require that organisms be reared to adulthood or that they be screened using elaborate testing procedures such as disease challenges or reproduction trials. With the advent of "marker-assisted" breeding, the selection process can be greatly accelerated by using genetic markers of observed phenotypic traits to select organisms at a young age without the need for phenotypic screening tests. This can increase the selection intensity and reduce the cost of a selective breeding program.

A prerequisite for marker-assisted breeding is identification of genetic markers for phenotypic traits of interest. The genetic marker need not be the specific gene responsible for a particular phenotypic trait, but it should be linked to it. In practice, this means that the marker is found together with the specific gene when the organism's DNA is fragmented. This implies that the marker is located in close proximity to the specific gene, or at least on the same chromosome.

The process of searching for markers involves the following steps:

- Fragmenting the DNA using restriction enzymes,
- Amplifying the DNA using PCR technology (PCR is a molecular biology technique that copies DNA fragments, increasing the initial number of fragments millions of times. A detailed explanation of PCR is available in the previous issue of this magazine, Vol. 2, Issue 4/5),
- Loading the DNA at the base of an electrophoretic gel,
- Applying an electric current to the gel to stimulate the negatively charged DNA fragments to migrate toward the positive pole,
- Staining the gel to identify the location of various DNA fragments, and
- Noting the bands created by faster migration of the smaller (lower molecular weight) fragments and slower migration of the larger fragments (Fig. 1).

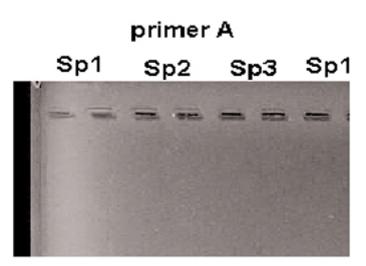


Fig. 1: Gel photography. Possible markers for these families are indicated in the yellow frames.

Usually many samples from different animals are run side-by-side on a gel. The vertical bands that form in the gel can be compared among samples. Bands that are common for all samples are called "monomorphic" and bands that are present in some animals but do not exist in others are called "polymorphic." Polymorphic bands are potential "markers."

There are many molecular biology techniques useful for finding these markers. The following is a presentation of three methods.

# **Random amplified polymorphic DNA (RAPD)**

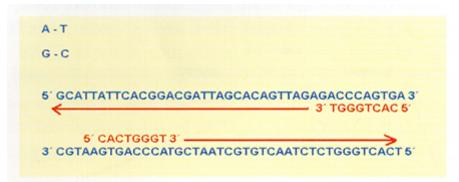


Fig. 2: RAPD uses a single primer (red), this primer bind to complementary sequences in DNA, and PCR produces a fragment amplification. The sizes of these fragments depend of the distance between to primers fixed in opposite strands.

This techniques utilizes a special Polymerase Chain Reaction (PCR) named AP-PCR (arbitrary Primed). RAPD can be run without pre-existing DNA sequence information. Only a single arbitrary primer is used, and the amplified fragments are genoma segments flanked by two complementary primer-binding sites in inverted orientation (Fig. 2). If the binding primer site has a modification in the sequence (mutation, deletion, insertion, etc.), the primer does not bind in this place and will produce a variation in the size of the band.

A set of primers is screened for reproducible polymorphisms. Different primers may identify different polymorphisms and have different reproducibility. Polymorphic bands may be potential genetic markers.

# Amplified fragment length polymorphism (AFLP)


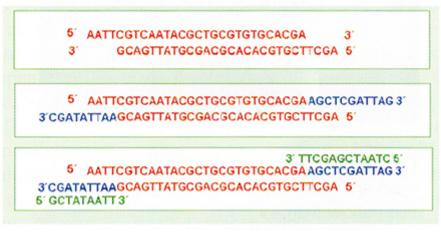



Fig. 3: Double strand of DNA with "sticking" ends, generated by enzyme digestion (red). Adapters binds to complementary sticking ends (blue). Primers (green) used in PCR reaction, fix to complementary sequence of adapters, and generate the amplification of this fragment.

The genomic DNA is digested with two restriction enzymes (restriction enzyme cut double-stranded DNA into specific sequences characteristic to each enzyme). One enzyme can cut a DNA sequence in a few select points, while another enzyme can cut the sequences in more places. Adapters (made of synthetic DNA) are ligated to the "sticky ends" of the genomic DNA, complementary to those of the restriction site (Fig. 3).

This results in the formation of a "hybrid" molecule, adapter-DNA-adapter. Since the sequences of the end fragments are known, we can elaborate primers that anneal (hybridize) to the ends' fragments. This DNA segment is then used as a template for the PCR reaction. The PCR products are separated on polyacrylamide gels, and the study of the bands is similar to that described in the RAPD technique. This technique, however, produces more bands than the RAPD and its results are more reproducible.

# Microsatellites

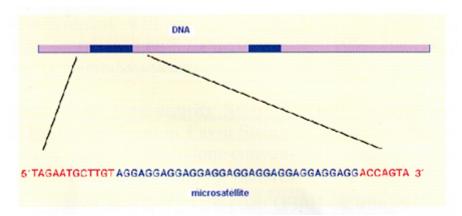

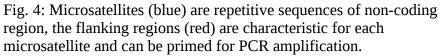

Microsatellites are repeated, non-coding sequences of DNA (do not correspond to genes), and are of variable length. Microsatellite variation can be detected by PCR using primers flanking the repeated DNA sequences (Fig.

#### 6/9/2024

4). There are two ways to identify microsatellite loci. One method is by searching through DNA sequence databases for sequences containing simple repeats.

Primers are then designed directly from the sequence data. The second approach is to digest the entire DNA and clone the fragments into a vector. The correct clone is identified using a probe. The DNA fragment is re-isolated from the clone in order to sequence the fragment. With luck, all of the microsatellite sequences with its flanking regions are determined. The primers are designed and the DNA extraction is submitted to the PCR, and then visualized electrophoretically.

Once the microsatellite markers are established, the utility of the marker is high. Fig. 5 shows the variation in the microsatellite sequence and its development in a gel. Different individuals within a given population have the same microsatellite (same sequence repeat), but the length of the sequence is different among these individuals. All posses the same flanking sequences and the specific primers will bind these flanking sites, but the fragment length obtained in the PCR will be different. It is these differences in length (and weight) which generate characteristic patterns of each individual.

Fig. 5: (a) Utilization of the same primers set (green arrows) in different individuals. All the individuals have the same microsatellite (represented by red and blue letters) flanked by the same sequence, but the length is different in each animal and may be different in homologous chromosomes of the same individual. (b) Schematic representation of a gel from chart (a).

From this brief description, it should be apparent that many of these new molecular biology technologies are not for everyone due to their complexity and cost. When embarking on a genetic selection program, it is strongly recommended to consider one or several of these tools. This, however, requires a commitment of both human and material resources in order to make the effort worth while.

(*Editor's Note: This article was originally published in the December 1999 print edition of the* Global Aquaculture Advocate.)

# Now that you've finished reading the article ...

... we hope you'll consider supporting our mission to document the evolution of the global aquaculture industry and share our vast network of contributors' expansive knowledge every week.

By becoming a Global Seafood Alliance member, you're ensuring that all of the pre-competitive work we do through member benefits, resources and events can continue. Individual membership costs just \$50 a year.

#### Not a GSA member? Join us.

#### Support GSA and Become a Member

#### Authors

• 💦 Sam Stern

Sam Stern

Desarrollo Industrial Bioacuatico (S.A.) A subsidiary of Conti Group Co. Ecuador

[99, 101, 46, 109, 111, 99, 46, 97, 115, 98, 105, 100, 64, 110, 114, 101, 116, 115, 115]

• 📄 Leonardo Galli

Leonardo Galli

Desarrollo Industrial Bioacuatico (S.A.) A subsidiary of Conti Group Co. Ecuador

## Share

- Share via Email
- **Share on Twitter**
- **f** Share on Facebook
- **in** Share on LinkedIn

## **Tagged With**

genetics selective breeding molecular biology Leonardo Galli Microsatellites Sam Stern

**Related Posts** 

Health & Welfare

## Antimicrobial resistance: complex, unavoidable

Resistance to antibiotics is unavoidable because it is an aspect of bacterial evolution. Resistant microflora have been observed in water that has not had recent contact with antimicrobial agents.

## Intelligence

## As ocean temperatures rise, so too will vibrio outbreaks

A study using a half-century of data has linked climate change and warming sea temperatures with an increase in illnesses from the common vibrio bacteria. Shellfish growers, fighting a particularly virulent strain of Vibrio parahaemolyticus, are changing their harvest protocols.

Health & Welfare

## **Bioinformatics: The ultimate breeding tool**

Bioinformatics, a true hybrid between molecular and computer sciences, changes the rules of both breeding and medicine with a novel approach.

Health & Welfare


## **Microsatellites: Genetic markers of choice**

Use of microsatellites is one of many marker approaches, but has proven to be the tool of choice across all commercial agriculture sectors.

## About The Advocate

The Responsible Seafood Advocate supports the Global Seafood Alliance's (GSA) mission to advance responsible seafood practices through education, advocacy and third-party assurances.

## Learn More

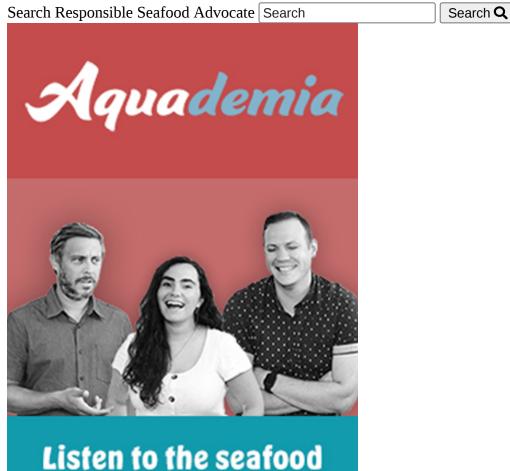

# Listen to the seafood industry's top podcast

## **Advertising Opportunities**

## 2022 Media & Events Kit

## Categories

Aquafeeds > Health & Welfare Health & Welfare > From Our Sponsors > Innovation & Investment > Intelligence > Responsibility > Fisheries > Artículos en Español >

## Don't Miss an Article

## Featured

- <u>Health & Welfare An update on vibriosis, the major bacterial disease shrimp farmers face</u>
- Intelligence A seat at the table: Fed By Blue team says aquaculture needs a stronger voice
- <u>Responsibility Quantifying habitat provisioning at macroalgae cultivation locations</u>

## **Popular Tags**

All Tags

#### Recent

- Fisheries Second Test: Another filler for the fisheries category
- <u>Fisheries Test: This is filler for the fisheries Category</u>
- Aquafeeds Test Article
- <u>Responsibility Study: Climate change will shuffle marine ecosystems in unexpected ways as ocean</u> <u>temperature warms</u>
- Health & Welfare Indian shrimp researchers earn a patent for WSSV diagnostic tool



# Listen to the seafood industry's top podcast

- <u>About</u>
- <u>Membership</u>
- <u>Resources</u>
- Best Aquaculture Practices (BAP)
- Best Seafood Practices (BSP)
- GOAL Events
- Advocate Magazine
- Aquademia Podcast
- <u>Blog</u>
- <u>Contact</u>

## Stay up to date with GSA

- 0
- f
- 🗙
- in
- 🕨

Copyright © 2024 Global Seafood Alliance All rights reserved. <u>Privacy</u> <u>Terms of Use</u> <u>Glossary</u>